
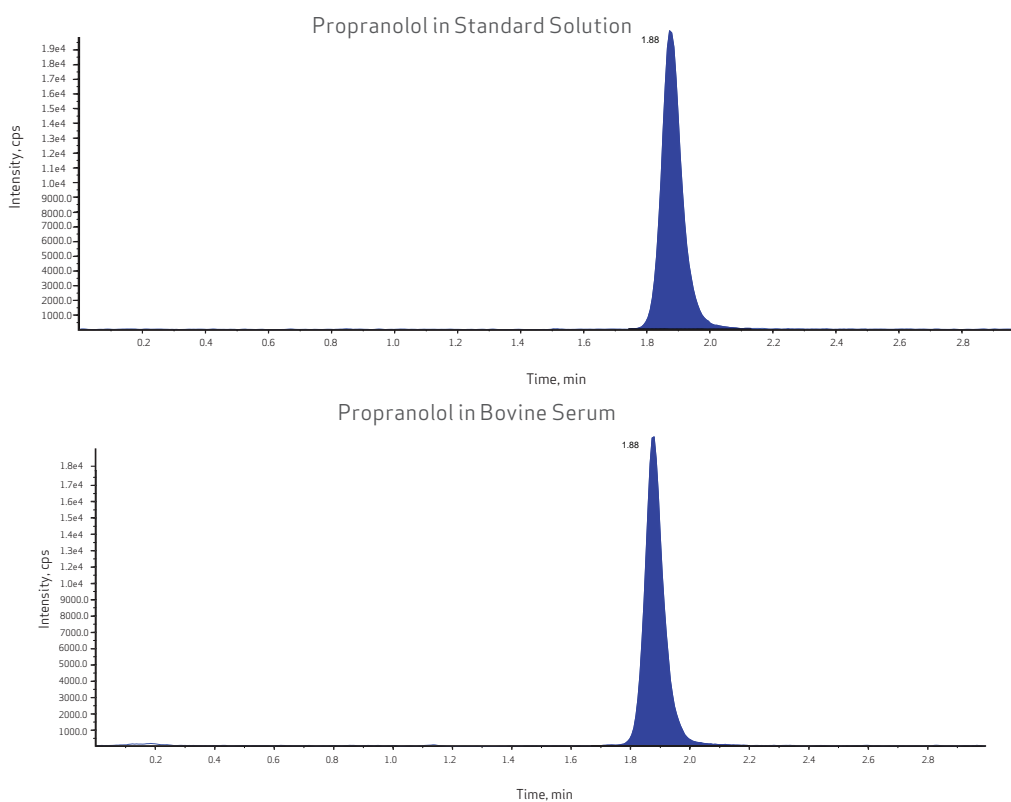


SiliaChrom® SB C18

High Performance Liquid Chromatography - Tandem Mass Spectrometry for Propranolol in Bovine Serum

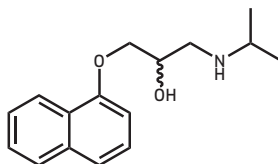

Conditions

Mass System: API 4000 Q TRAP tandem mass spectrometer

Mass Condition: Multiple reaction monitoring (MRM) mode; electrospray ionization (ESI) in the positive ion mode; m/z 409.1→238.1 was used to quantify amlodipine

Column: SiliaChrom SB C18, 5 µm

Column Size: 2.1 x 150 mm

Mobile Phase: 0.01M ammonium acetate solution containing 0.1% formic acid/acetonitrile (40/60, v/v)

Flow Rate: 0.5 mL/min

Temperature: 25°C

Injection volume: 5 µL

© Copyright 2010 SiliCycle Inc. All rights reserved. Subject to change without notice.